


Socket Repair Membrane


Technique For Facial Plate Defects

Atraumatically Extract Tooth Using A Flapless Technique

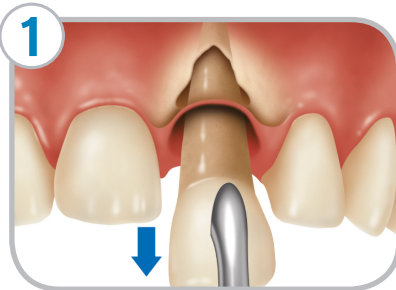

1 Prepare the bone graft material (e.g. Puros® Particulate) according to the instructions for use. Utilizing an atraumatic, flapless technique, carefully extract the tooth and completely debride the socket.


Shape & Prepare Membrane

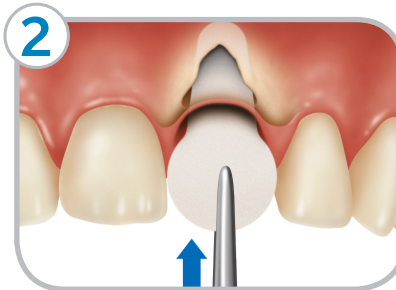

2 Remove the Zimmer® Socket Repair Membrane from the packaging. Insert the small end of the pre-cut membrane into the socket, ensuring that it extends over the facial defect both laterally and apically. Establish how much of the membrane's large end will be needed to fold closely

over the top of the socket to facilitate closure after grafting. Remove the membrane and trim as needed.


Line Facial Defect With Membrane

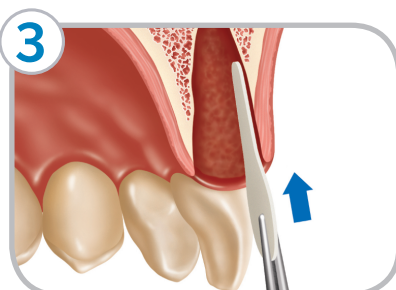

3 Insert the small end of the membrane into the socket and against the facial tissue.

Caution: To help maintain blood supply to the remaining facial plate, do not lift the periosteum off the bone.


Augment Socket With Bone Graft Material

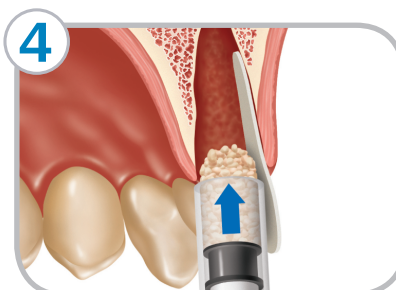

4 Augment the socket with graft material using a sterile syringe or sterile applicator. To fill the syringe, pull back on the plunger and gently press its tip into the hydrated graft material. Repeat until the plunger is filled. After placing bone graft material into socket,


firmly compress with a sterile instrument. This will prevent voids at the apex of the socket and push the facial tissue labially for better ridge shape.

Close Socket with Membrane

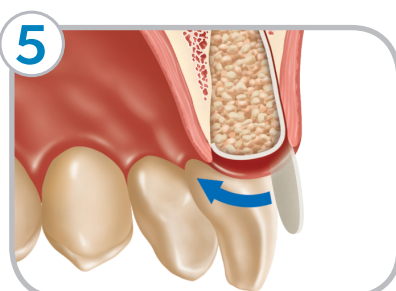

5 Gently close the socket by folding the membrane over the bone graft material until the membrane approximates with the palatal/lingual tissue. The membrane should be stable due to the graft material pressing it against the facial plate.

Suture Membrane to Soft Tissue

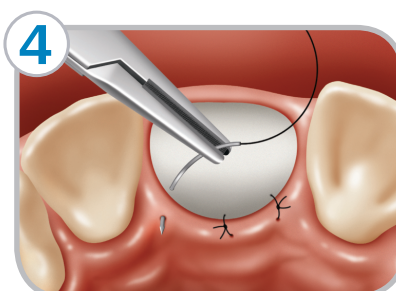

Place at least 2 to 3 interrupted 5-0 absorbable sutures to secure the membrane to the palatal/lingual tissue.

Note: Typically, no sutures on the labial/buccal aspect are needed since the membrane will be held by pressure of the graft material.

SOCKET REPAIR MEMBRANE TECHNIQUE FOR FACIAL PLATE DEFECTS

The Zimmer Socket Repair Membrane is designed to assist wound healing in alveolar facial plate repair and residual ridge preservation following atraumatic, flapless single-root tooth extraction.

- Socket grafting can help to preserve bone volume for implant placement and improve aesthetic results.¹
- The socket repair procedure is a flapless technique designed to preserve natural soft tissue architecture and vascularity.²
- Guided bone regeneration is performed inside the socket by lining the defect with the bioabsorbable collagen membrane prior to filling with bonegraft material.² The membrane is designed to stabilize the graft material and block fibrous tissue ingrowth.²
- An external extension of the membrane is folded over the top of the socket and sutured to the palatal/lingual soft tissue to help contain the graft material² and aid in stabilizing the clot for healing.³
- As with any socket grafting procedure, minor loss of coronal graft particles may normally occur until soft tissue healing provides complete coverage.²
- Membrane is usually completely resorbed 26-38 weeks following surgery.*

Suggested surgical components for this procedure:

- | | |
|--|--|
| <ol style="list-style-type: none"> 1 Zimmer Socket Repair Membrane 2 Bone graft material from Zimmer Biomet Dental (e.g. Puros Particulate) 3 1 cc open-bore syringe or sterile applicator to deliver the bone graft material | <ol style="list-style-type: none"> 4 Tissue forceps 5 Scissors to trim membrane 6 5-0 absorbable suture |
|--|--|

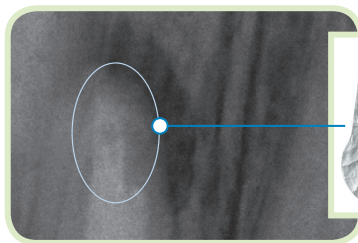
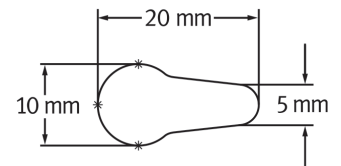


Figure 1 Periapical pathology caused loss of buccal plate of bone.

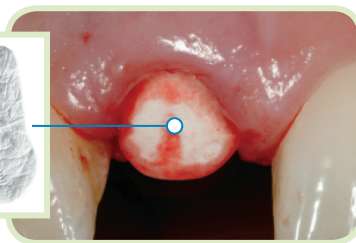


Figure 2 Membrane placed in socket.

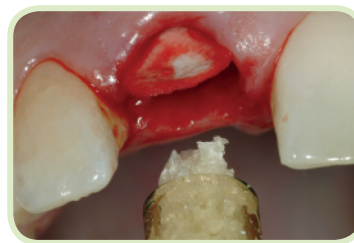


Figure 3 Puros Particulate delivered into socket.

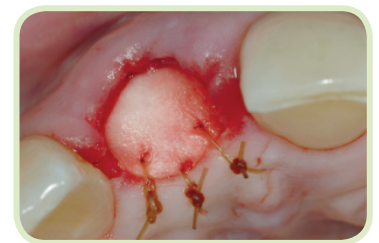


Figure 4 Membrane sutured in place.



Figure 5 Tissue healing at two weeks.



Figure 6 Final aesthetic result.

Clinical photos ©2010 Dennis Tarnow, DDS, and Richard Smith, DDS. All rights reserved. Individual results may vary.

*When not exposed the resorption rate is 26-38 weeks; if left exposed, resorption time is shorter.

1. Sclar AG. Ridge preservation for optimum esthetics and function: the "Bio-Col" technique. *Postgraduate Dentistry*. 1999;6:3-11.
2. Elian N, Cho SC, Froum S, Smith RB, Tarnow DP. A simplified socket classification and repair technique. *Pract Proced Aesthet Dent*. 2007;19:99-104, quiz 106.
3. Bunyaratavej P, Wang HL. Collagen membranes: a review. *J Periodontol*. 2001; 72:215-229.

Ordering Information: Zimmer Socket Repair Membrane • Catalog Number 0154

Contact us at 1-800-342-5454 or visit zimmerbiometdental.com

Zimmer Biomet Dental
Global Headquarters
4555 Riverside Drive
Palm Beach Gardens, FL 33410
Phone: +1-561-776-6700
Fax: +1-561-776-1272

Unless otherwise indicated, as referenced herein, all trademarks are the property of Zimmer Biomet; and all products are manufactured by one or more of the dental subsidiaries of Zimmer Biomet Holdings, Inc., and distributed and marketed by Zimmer Biomet Dental (and, in the case of distribution and marketing, its authorized marketing partners). Zimmer Socket Repair Membrane is manufactured by Collagen Matrix, Inc. Puros products are manufactured by RTI Surgical, Inc. For additional product information, please refer to the individual product labeling or instructions for use. Product clearance and availability may be limited to certain countries/regions. This material is intended for clinicians only and does not comprise medical advice or recommendations. This material may not be copied or reprinted without the express written consent of Zimmer Biomet Dental. ZBINST0001 REV A 05/17 ©2017 Zimmer Biomet, All rights reserved.

